

SKB041-3

PATIENT TRANSPORTATION TROLLEY



Healthward Tuk

Healthier your life

SKB041-3











FEATURES AND HIGHLIGHTS

Technical parameters

X-ray Backrest

External size (LxW):	2120x840 mm
Height from floor to bed platform:	550-860 mm
Height from floor to lowered guardrail:	590-900 mm

X-ray Platform

External size (LxW):	2120x840 mm
Height from floor to bed platform:	600-910 mm
Height from floor to lowered guardrail:	640-950 mm

Back plate adjustment:	0-85° (±5°)
Knee-rest adjustment:	0-30°(±5°)
Reverse Trendelenburg:	0-15°
Trendelenburg:	0-15°

Technical configuration

☑ Usa hydraulic pump	2 pcs
✓ Crank system	1 set
☑ Back rest air spring	2 pcs
☑ 8" German castors,central locking	4 pcs
☑ Emergency push handle	1 set
☑ Directional wheel	1 pc
✓ Folding side rail	1 set
☑ Ф74mm wall bumper wheel	4 pcs
✓ Iv pole previsions	4 pcs
✓ Mattress	1 pc
☑ IV pole	1 set
Central locking break pedal	2 sets
☑ Directional wheel release pedal	2 pcs
☑ X-ray translucent bed platform	1 set







Metal material

Baosteel from Fortune Global 500 Panasonic Robotic ensure 360° full smooth wedlding.



Plastic

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.



Steels | Painting

11 Process epoxy painting, ASTM testing anti-baterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.



Brake

Break spring keep qualified after 10,000 times uninteruped break.



Hydraulic pumps

Height adjustment, trendelenburg and revised tendelenburg through two USA hydraulic pumps.



Center wheel

Can be engaged to facilitate direction control, enabling users to navigate the stretcher easily in hectic and dynamic environments.




Rising


Hydraulic treadle for rise.



Falling

① Leg section lows down to make head highter.

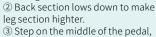

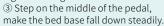



Castors

Four double side casters of 6 inch diameter, simultaneously lockable by the step pedal at both sides. TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts.



Mattress

Waterproof PU mattress, convenient for transfer patient.



Side rail

Foldable aluminium alloy guardrail



Bumper wheel

Bumper wheel detachable, max. protect the bed and patient during transportation



Drainage Hooks

Mobile plastic drainage hooks.



IV pole

Stainless Steel Column with plastic hook.

















Design

Research

Manufacturing

Maintenance

Service



USA

Saikang Medical Technology, Inc. 17800 Castleton St Ste 180, City of Industry, 91748, California. Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC Truda 174, 454080 Chelyabinsk. Tel.: +7 (351) 2145566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd. Viraj Complex Unit 3, Mombasa Road, P.O. Box 14262-00800, Nairobi. Kenya. Tel.: +254733704690

Middle East

Saikang Medical Supplies LLC X22-S03 England Cluster, International City, Dubai. United Arab Emirates. Tel.: +20 128 8962497

Healthward uk

Chengdu Saikang Medical Equipment Co., Ltd. STOCK CODE: 870098

- No.1, 5th Floor, Building 5, No.360, Wuke West 5th Road
- +86 18021231901 · \$\sum_2\$ 24/7
- export@health-ward.com
- www.health-ward.com







Company video

General catalogue